


Julian Bergmann^[1], Alexander Pikhtelev^[2], Samuel Ayet San Andres^[3], Timo Dickel^[1,3], Jens Ebert^[1,4], Hans Geissel^[1,3], Christine Hornung^[1], Johannes Lang^[1], Wayne Lippert^[1], Christian Lotze^[1], Wolfgang Plaß^[1,3], Ann-Kathrin Rink^[1] and Christoph Scheidenberger^[1,3]

^[1] JLU Gießen, Deutschland — ^[2] Institute of Energy Problems of Chemical Physics, Russian Academy of Sciences, Moscow, Russia — ^[3] GSI Darmstadt, Deutschland — ^[4] TU Darmstadt, Deutschland

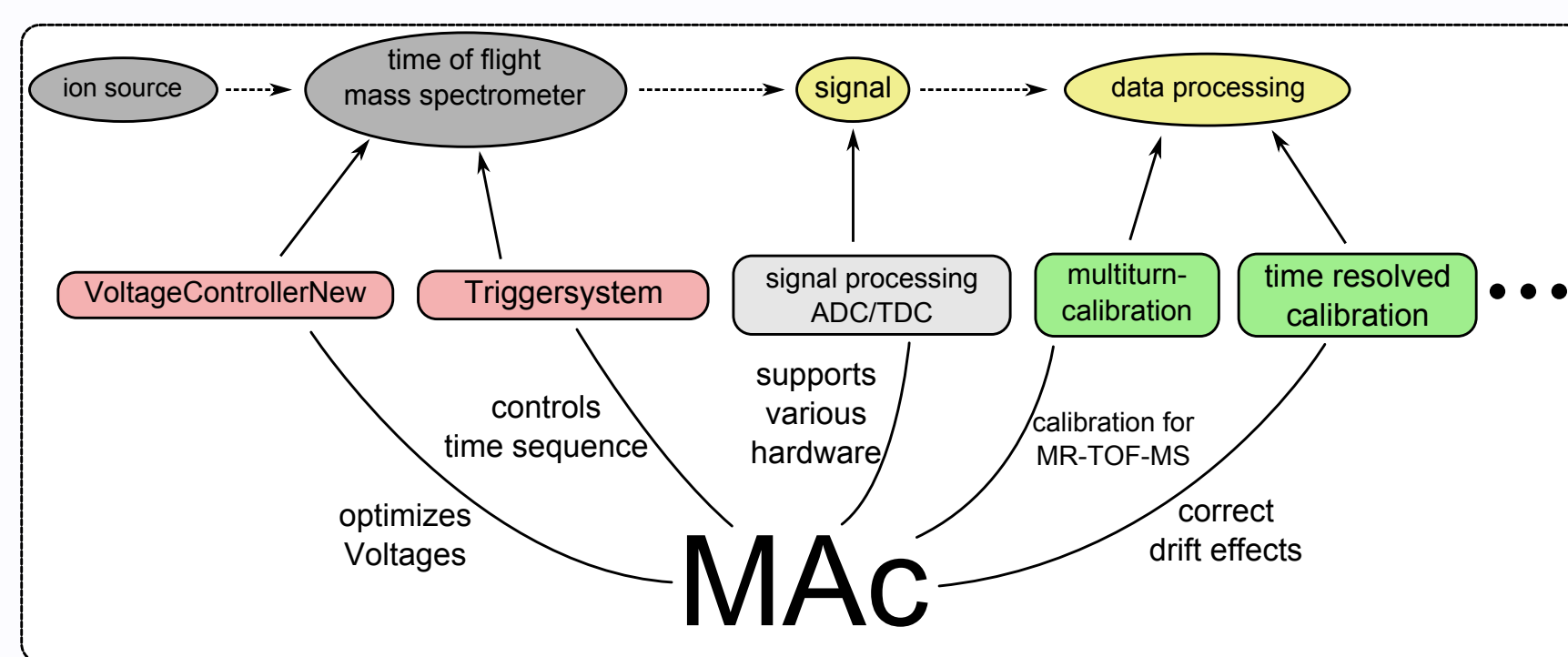

Features:

- + ADC/TDC for large dynamic range
- + time sequence controller
- + multi-turn calibration
- + voltage optimizer
- + time resolved calibration
- + spectra accumulation
- + automatic peak detection
- + import/export various file formats
- + data operations (e.g. smoothing)

2. Enlarging Dynamic Range Using TDC and ADC Simultaneously.

TDC:

- + detection of single ions
- dead time effect

ADC:


- + detection of high signal rates
- low intensity signals lost

Combined:

- + dynamic range > 4 orders of magnitude
- + detection of high and low signal rates

MAC:

- common user interface
- common calibration
- common file format and analysis tools

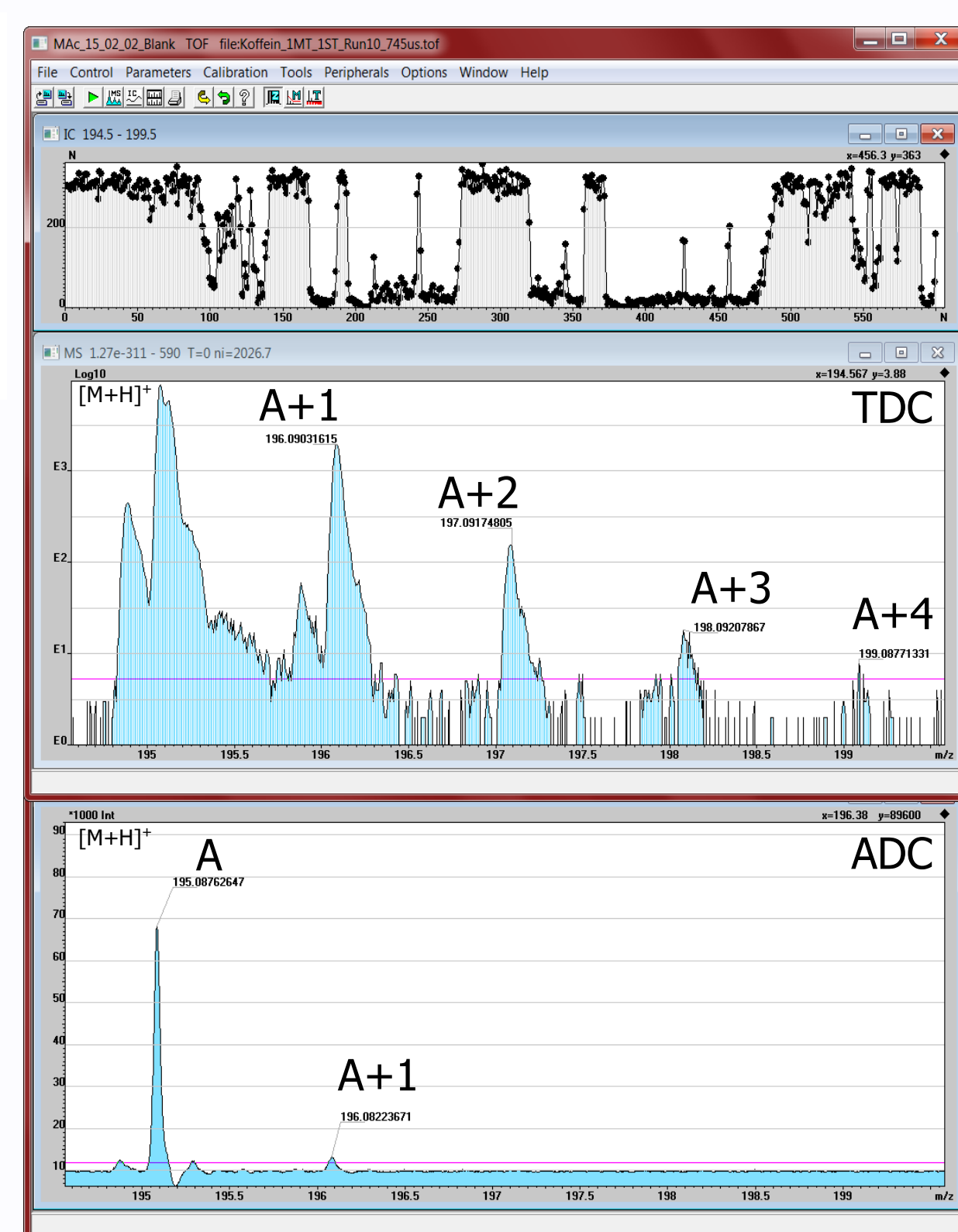

Example: Coffein (10^{-4} mol)

TDC (top):

- + A+1 to A+4

ADC (bottom):


- + A and A+1

Device	Isotope	MAc	Literature	Rel. Deviation
ADC	m195 / m196	9,374 79	9,457 79	0,008 78
	m196 / m197	11,7599	11,7515	$7,1480 \cdot 10^{-4}$
TDC	m197 / m198	5,629 95	19,1104	0,705 40
	m198 / m199	3,404 73	33,768	0,899 17


excellent isotope distribution
deviations from literature values
due to unresolved isobars

3. Time Sequence Controller

- controls time sequences for RF ion traps, analyzer electrodes etc.
- multiple modes to set delays and duty cycles

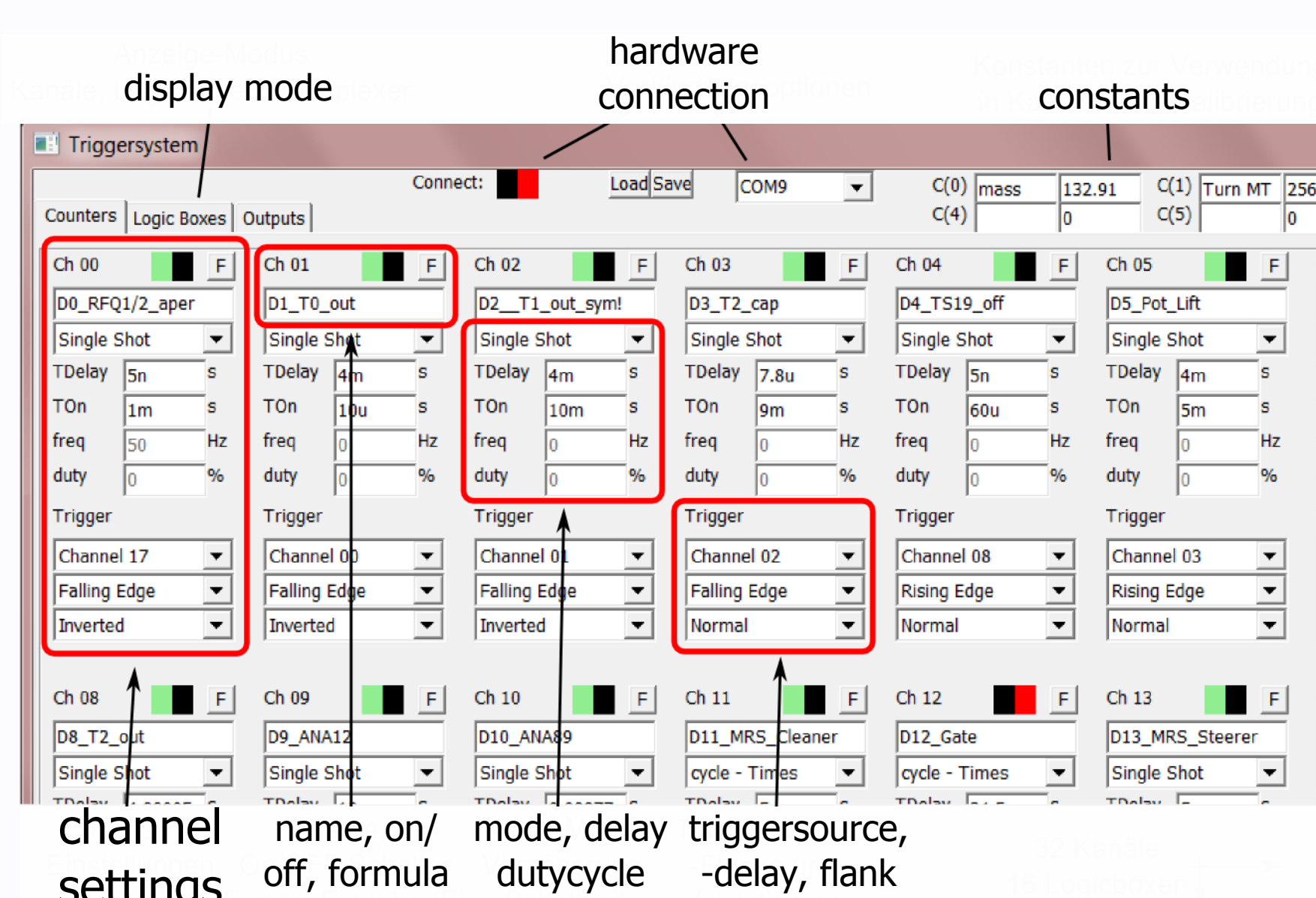

Example of automatic scaling

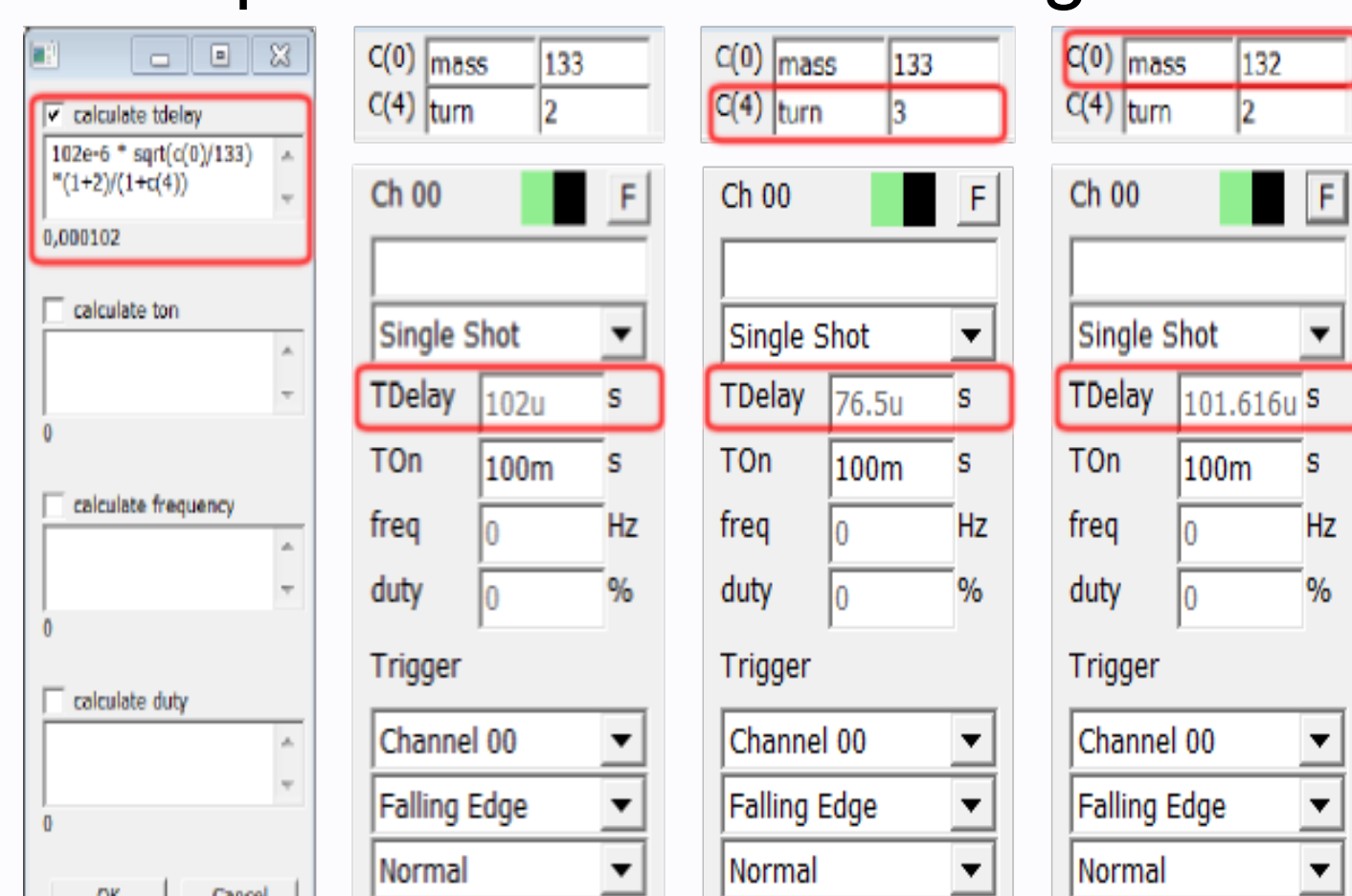

- + links triggers for 32 channels to each other or up to 16 logic channel combinations
- + uses formulas to calculate delay and duty cycles
- + link to other channel's times
- + use up to 8 constants

- + set times can be used in other MAC parts
- calibration

4. Multiturn Calibration

Calibration Formula

$$m = a \frac{(t - t_0)^2}{(1 + b \cdot N)^2}$$

$$a = \frac{2q \cdot U}{l_{detec}^2} \quad b = \frac{l_{turn}}{l_{detec}}$$


t_0 : electronic delay time

- + ≤3 calibrants → analytical solution
- + >3 calibrants → least square deviation minimizer
- + uses time sequence controller for trap extraction time
- + use and display for different turns in multiturn calibration
- + mass range greatly increased


Example: Xe & SF6

- + between 126 and 130 turns
- + precision of online identified masses of $4 \cdot 10^{-7}$

Time (µs)	Turns	Mass	MAc (u)	Mass	Lit (u)	Lit-MAc	Nuklid
49,3071469	130	127,9029259	127,9029822	4,40 · 10 ⁻⁷	128	128	¹²⁸ Xe ⁺
46,5327600	126	135,906564	135,9066657	7,48 · 10 ⁻⁷	136	136	¹³⁶ Xe ⁺

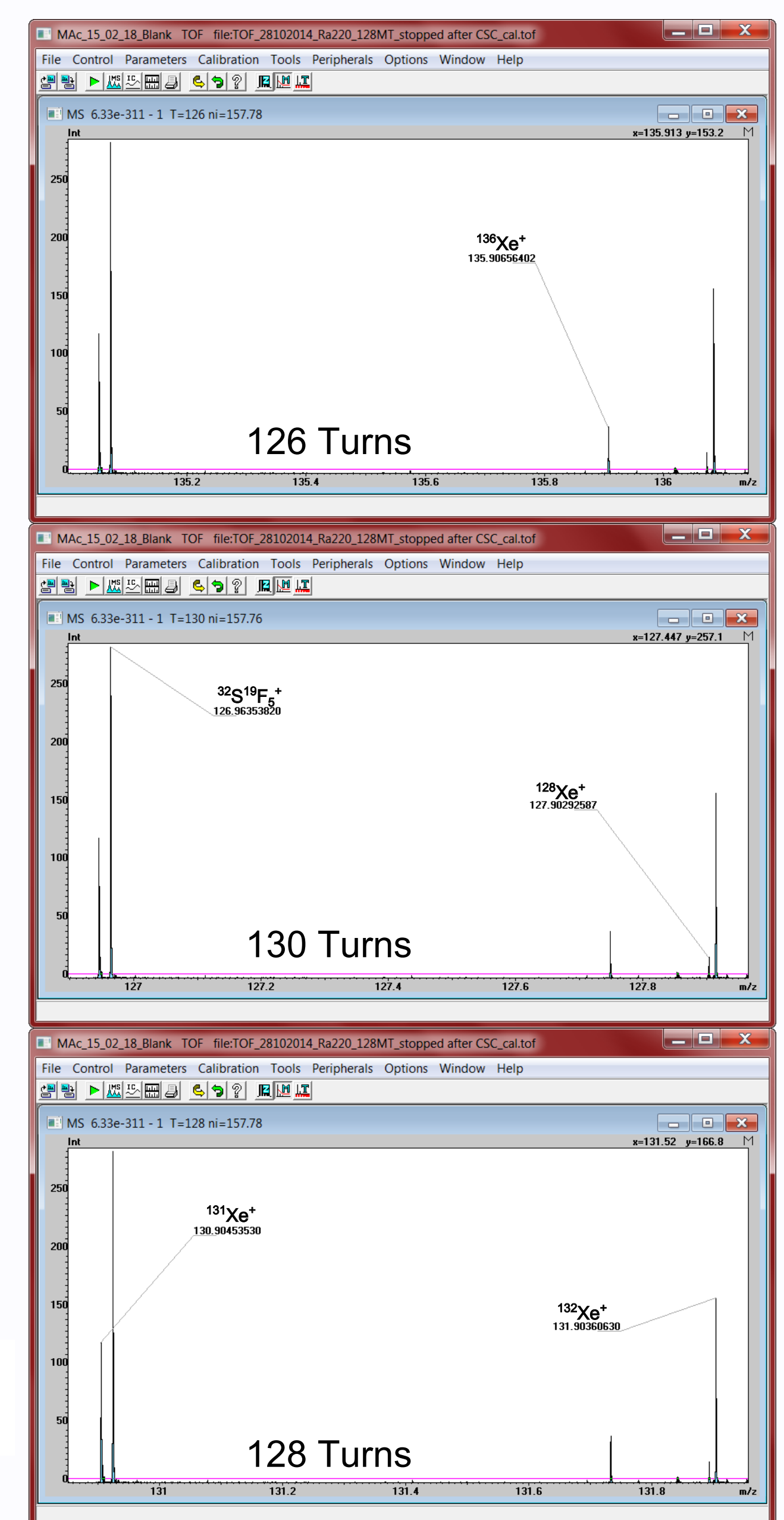

5. Voltage Optimizer

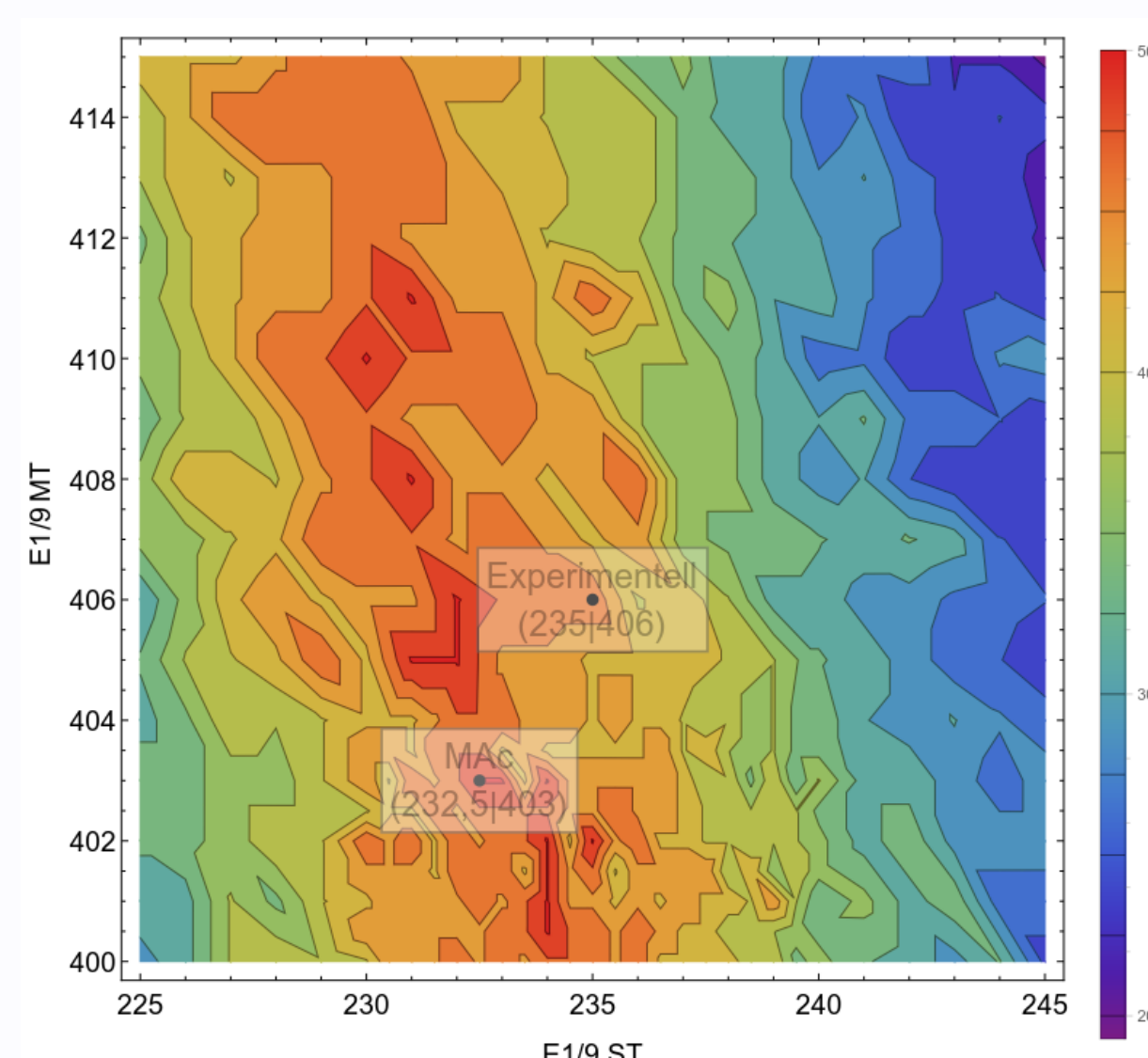

- connects to all available power supplies
- free choice for amount of channels
- from adjustable MIN to MAX all voltage combinations are scanned
- automatic measurement taken for each combination
- + up to 20 measurements per minute
- automatic optimization algorithms to be implemented

Conclusion:

Indispensable for preparing optimized off- and online measurements


6. Time resolved calibration

Problem:

- voltages can vary during measurement
- thus flight length varies during turns
- mass lines shift in time

Solution:

- + each spectra gets calibrated → calibrant stays fixed in mass
- + automatic detection of calibrant
- + correction of 2000 spectra in under 1 second
- + recalibrated data can be displayed and used inside MAC
- + ideally suited for long term and high resolution measurements

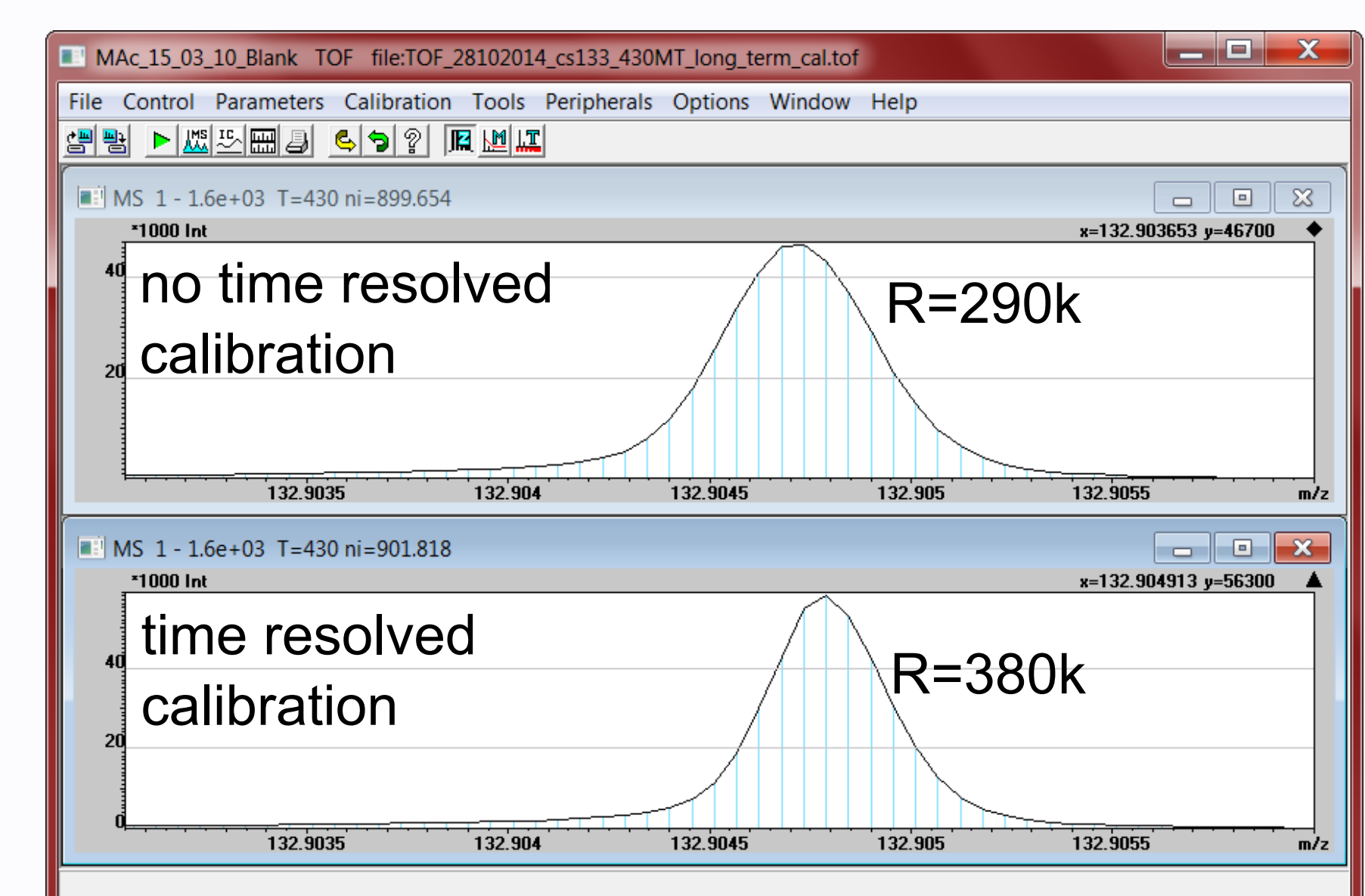

Poster Information

Presentation at DPG Heidelberg 2015
Poster MS 6.3
Discussion: March 24, 17:00 until 19:00